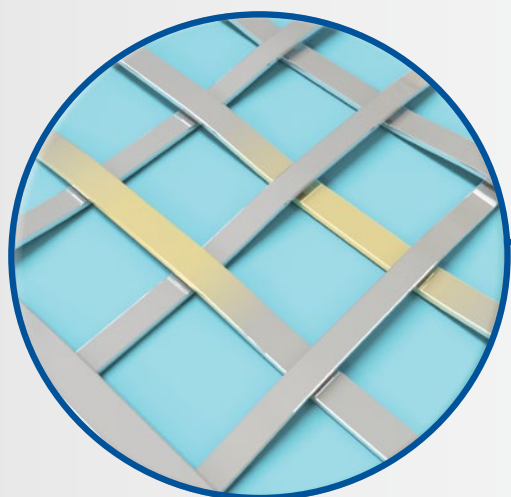


5F and 6F Now Available

GET MORE

Control
Back-up Support
Space to Work

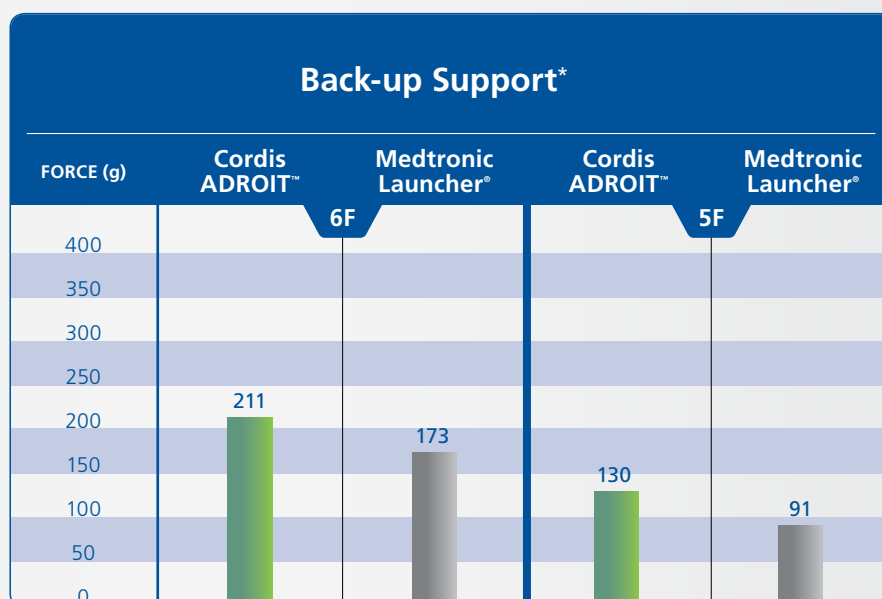


INNOVATIVE HYBRID DESIGN FEATURES TWO DIFFERENT SIZES OF WIRE BRAID

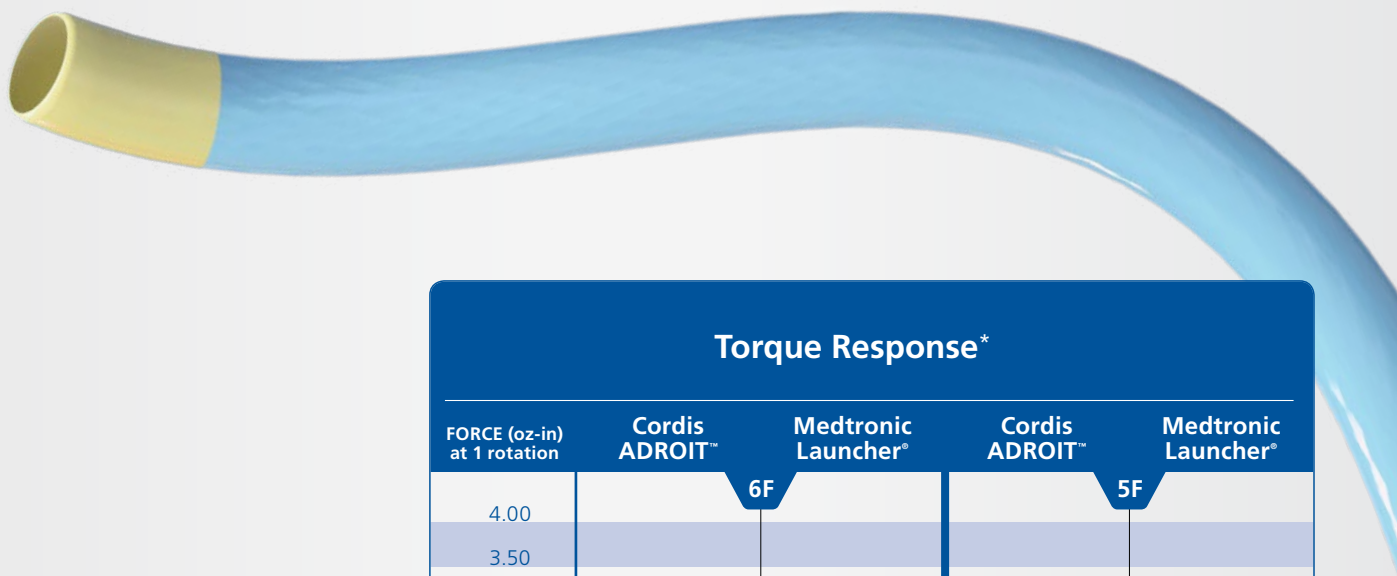
MORE BACK-UP SUPPORT

Engineered for strong back-up support.

- Innovative wire braiding with “flat-on-flat” and “2-over-2-under” design—providing design flexibility, kink resistance, and back-up support

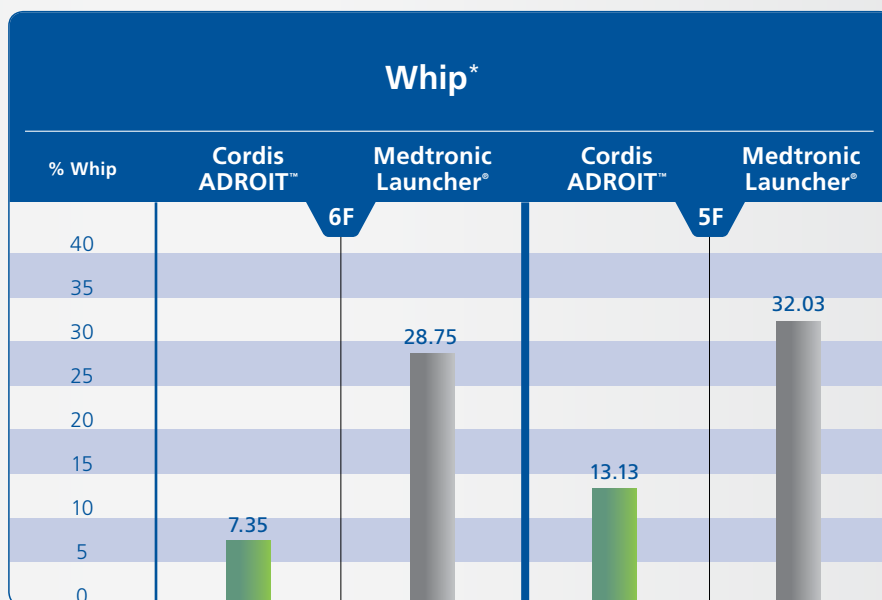
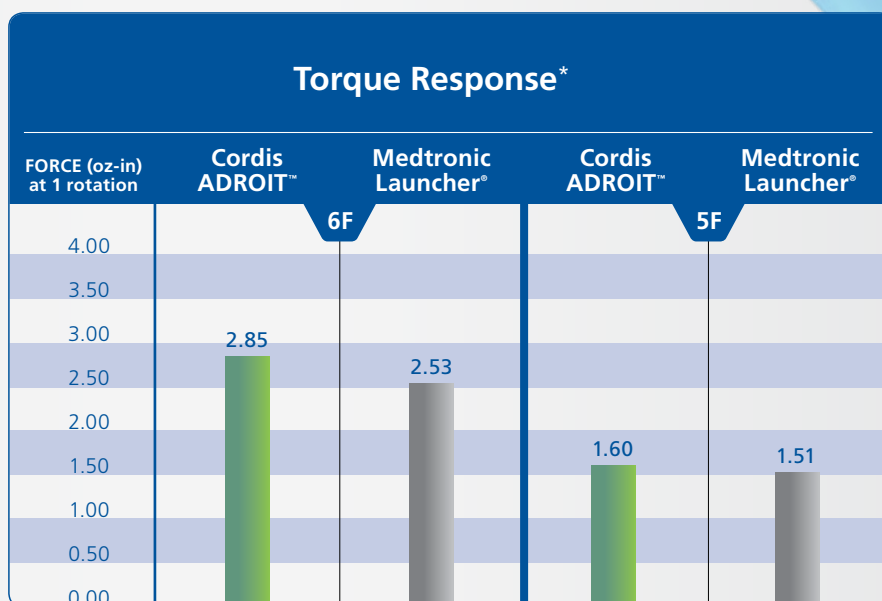


MORE CONTROL



Gain greater positioning confidence through:

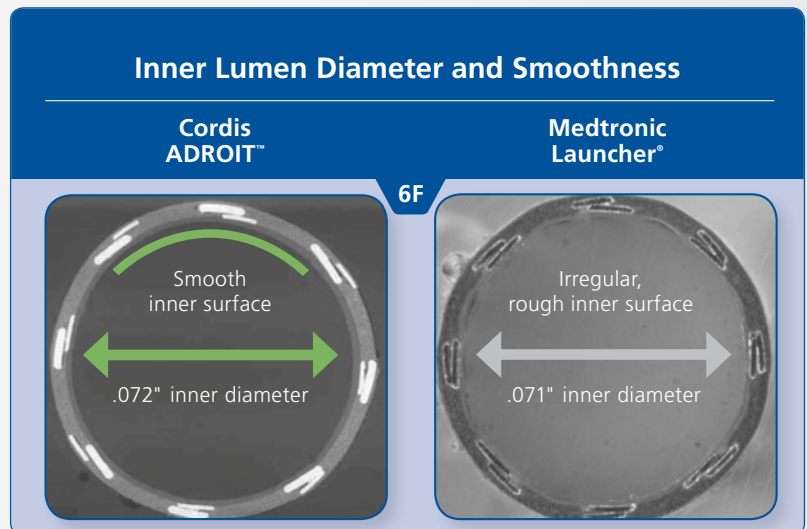
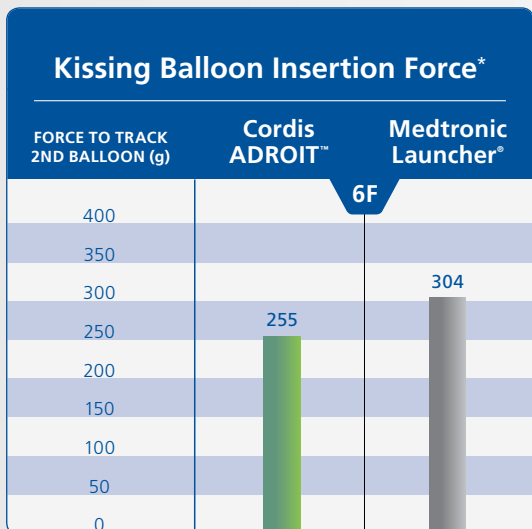
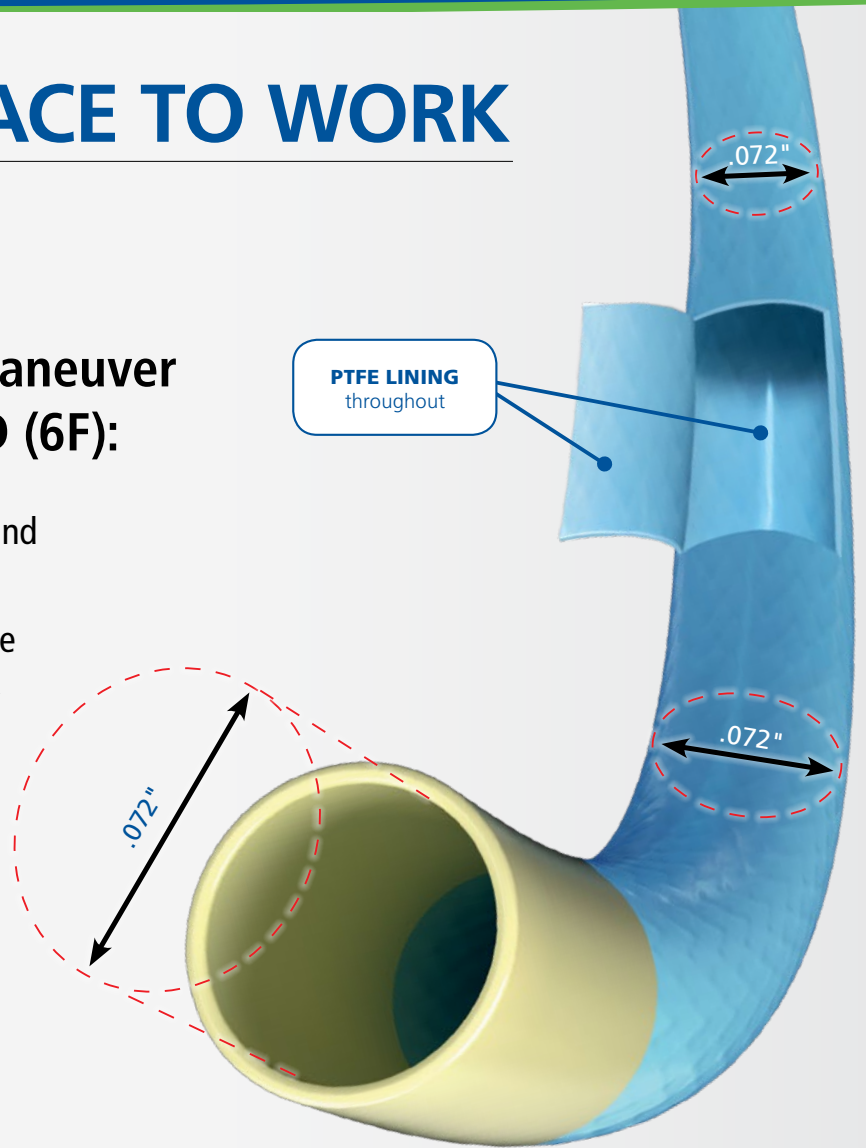
- Hybrid braided wire technology that delivers excellent torque response
- An advanced multi-segment design that improves kink resistance and minimizes whip



MORE SPACE TO WORK

Experience more room to maneuver with the largest available ID (6F):

- .072" ID to improve device compatibility and provide better visualization†
- Easier to perform kissing balloon technique
- Innovative hybrid braided wire technology enables larger lumen with optimal back-up support
- PTFE lining provides lubricious lumen for smoother delivery
- TRUELUMEN® Technology ensures consistent ID from hub to tip, for added confidence



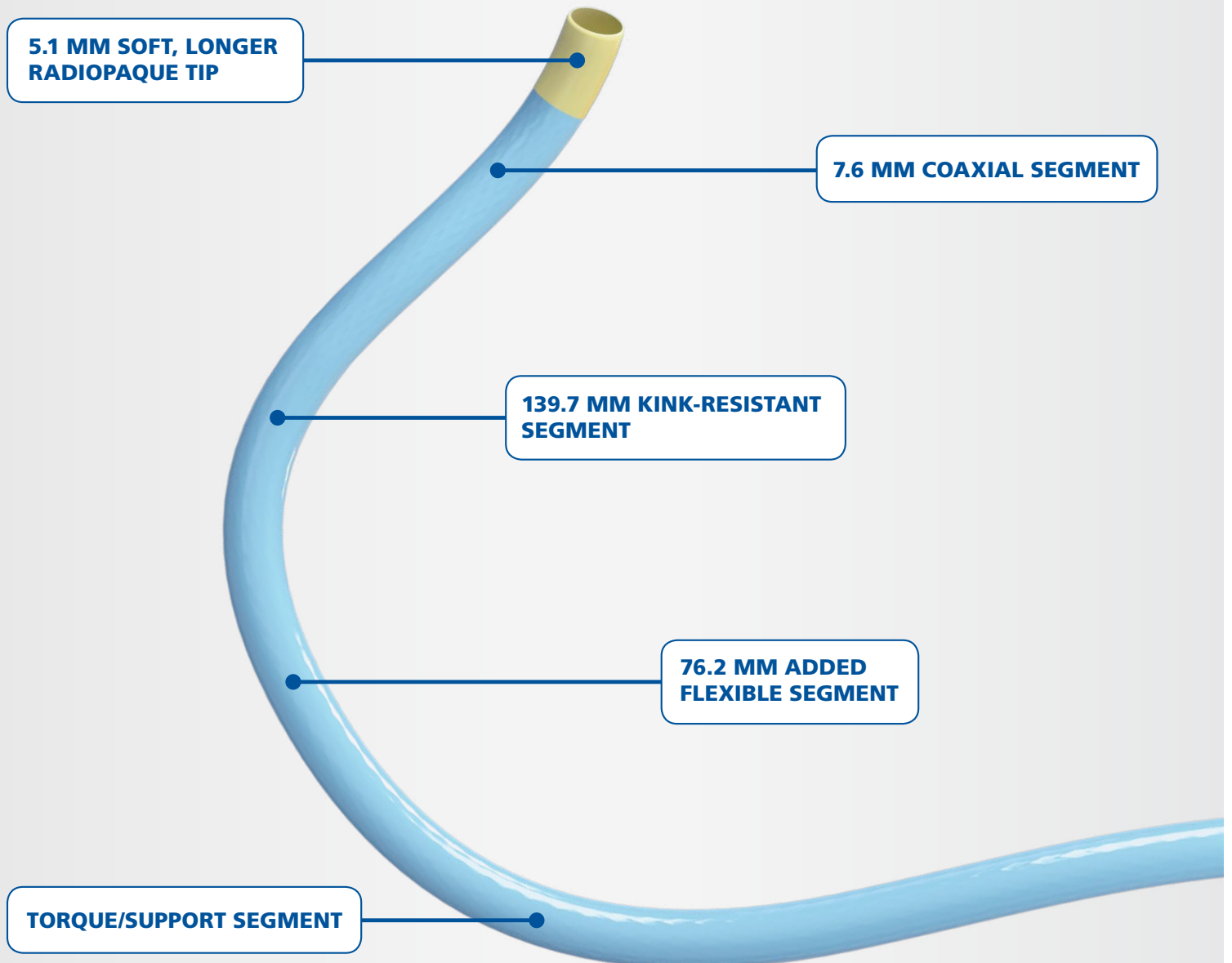
Innovative hybrid braided wire technology enables larger lumen diameter (.072") than Medtronic Launcher.®

†Cordis ADROIT™ Guiding Catheter 5F ID=.058"

DESIGNED FOR ATRAUMATIC PLACEMENT

Engineered to avoid vessel wall damage via:

- Longer, soft radiopaque tip allows active cannulation and provides enhanced visibility
- An advanced multi-segment construction improves flexibility and minimizes whip



GET MORE



Control

Gain more control over your procedure through improved kink resistance, decreased whip percentage, and excellent torque response.

Back-up Support

Feel more back-up support through the innovative hybrid braided wire technology and a wide range of extra back-up shapes.

Space to Work

Experience more space to work with improved device compatibility and better visualization using the largest ID available.

ADROIT™ Guiding Catheter from Cordis

Speak with your Cordis representative about getting more out of your guiding catheter.

Cordis
PART OF THE *Johnson & Johnson* FAMILY OF COMPANIES

*Limited internal testing data on file. Test Report 100323183, 100061993. Cordis Corporation. Tests were not powered to show statistical differences.

The trademarks used herein are trademarks of their respective owners.

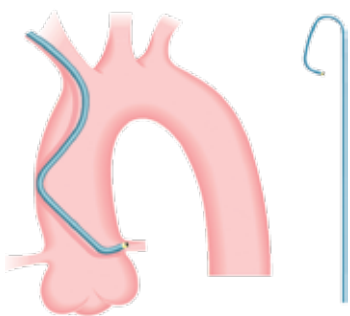
For EMEA healthcare professionals only. As part of the Cordis policy of continuous product development, we reserve the right to change product specifications without prior notification.

Important information: Prior to use, refer to the Instructions for Use supplied with this device for indications, contraindications, side effects, suggested procedure, warnings, and precautions.

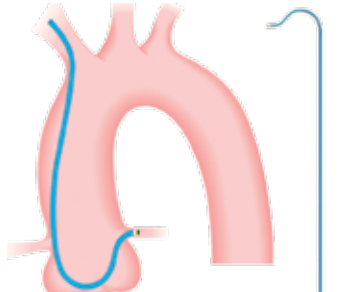
© Johnson & Johnson Medical NV/SA 2015 – EU737 5/15

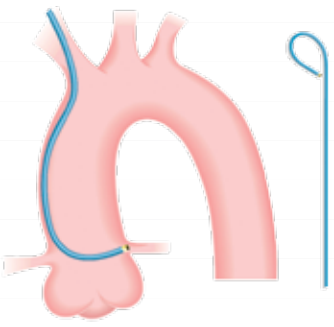
Cordis **ADROIT**[™] Guiding Catheter

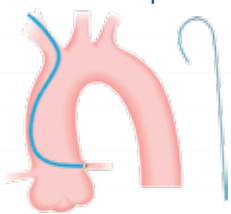
A large portfolio of shapes and configurations are available to meet varying procedural and anatomical demands.

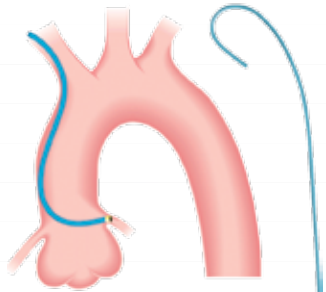
Left Coronary Artery	Shape	Description	PRODUCT CODE	
			6F	5F
 <p>Judkins Left</p>	JL 3		672-000-00	558-000-00
	JL 3.5		672-002-00	558-002-00
	JL 3.5 SH	2 Side Holes	672-003-00	
	JL 3.5 ST	Short Tip	672-012-00	
	JL 4		672-004-00	558-004-00
	JL 4 (90 cm)		672-004-90	
	JL 4 (125 cm)		672-040-25	
	JL 4 SH	2 Side Holes	672-005-00	
	JL 4 ST	Short Tip	672-014-00	
	JL 4 LBT	Long BRITE TIP [®]	672-004-0L	
	JL 4	ECOPAK (5 units)	672-004-0E	
	JL 4.5		672-006-00	558-006-00
	JL 4.5 ST	Short Tip	672-016-00	
	JL 5		672-008-00	558-008-00
	JL 6		672-010-00	

Left Coronary Artery	Shape	Description	PRODUCT CODE	
			6F	5F
 <p>Judkins Curved Left</p>	JCL 3.5		672-026-00	558-026-00
	JCL 3.5 SH	2 Side Holes	672-027-00	
	JCL 4		672-028-00	

Left Coronary Artery	Shape	Description	PRODUCT CODE	
			6F	5F
 <p>Amplatz Left</p>	AL .75		672-034-00	558-034-00
	AL .75 SH	2 Side Holes	672-035-00	
	AL 1		672-036-00	558-036-00
	AL 1 SH	2 Side Holes	672-037-00	
	AL 1 ST	Short Tip	672-044-00	
	AL 1.5		672-038-00	558-038-00
	AL 2		672-040-00	558-040-00
	AL 2 SH	2 Side Holes	672-041-00	
	AL 2 LBT	Long BRITE TIP [®]	672-040-0L	
	AL 3		672-042-00	558-042-00

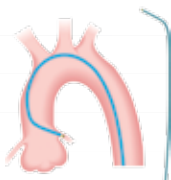
Left Coronary Artery (cont.)	Shape	Description	PRODUCT CODE	
			6F	5F
 <p>Extra Back-up</p>	XB 2		672-464-00	
	XB 2.5		672-465-00	
	XB 2.75		672-730-00	
	XB 3		672-052-00	558-052-00
	XB 3 SH	2 Side Holes	672-053-00	
	XB 3.25		672-400-00	558-400-00
	XB 3.5		672-054-00	558-054-00
	XB 3.5 (90 cm)		672-054-90	
	XB 3.5 (125 cm)		672-054-25	
	XB 3.5 SH	2 Side Holes	672-055-00	
	XB 3.5	ECOPAK (5 units)	672-054-0E	
	XB 3.5 LBT	Long BRITE TIP®	672-054-0L	
	XB 4		672-056-00	558-056-00
	XB 4 SH	2 Side Holes	672-057-00	
	XB 4.5		672-058-00	558-058-00
	XB 4.5 SH	2 Side Holes	672-059-00	

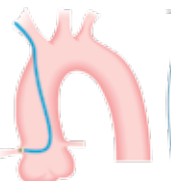
 <p>Extra Back-up LAD</p>	XBLAD 3		672-066-00	
	XBLAD 3 SH	2 Side Holes	672-067-00	
	XBLAD 3.5		672-060-00	558-060-00
	XBLAD 3.5 SH	2 Side Holes	672-061-00	
	XBLAD 4		672-062-00	558-062-00
	XBLAD 4 SH	2 Side Holes	672-063-00	
	XBLAD 4.5		672-064-00	


 <p>Extra Back-up Contralateral</p>	XBC 2.5		672-732-00	
	XBC 2.75		672-734-00	
	XBC 3		672-070-00	558-070-00
	XBC 3 SH	2 Side Holes	672-071-00	
	XBC 3.25		672-736-00	
	XBC 3.5		672-072-00	558-072-00
	XBC 3.5 SH	2 Side Holes	672-073-00	
	XBC 4		672-074-00	558-074-00
	XBC 4 SH	2 Side Holes	672-075-00	
	XBC 4.5		672-076-00	
	XBC 4.5 SH	2 Side Holes	672-077-00	

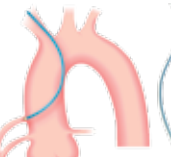
Right Coronary Artery	Shape	Description	PRODUCT CODE	
			6F	5F
 <p>Judkins Right</p>	JR 3.5		672-080-00	558-080-00
	JR 3.5 SH	2 Side Holes	672-081-00	
	JR 4		672-082-00	558-082-00
	JR 4 (90 cm)		672-082-90	
	JR 4 (125 cm)		672-082-25	
	JR 4 SH	2 Side Holes	672-083-00	
	JR 4 ST	Short Tip	672-090-00	
	JR 4 LBT	Long BRITE TIP®	672-082-0L	
	JR 4	ECOPAK (5 units)	672-082-0E	
	JR 5		672-084-00	558-084-00

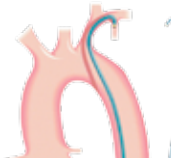
Right Coronary Artery (cont.)		Shape	Description	PRODUCT CODE	
				6F	5F
Judkins Curved Right		JCR 4		672-098-00	
Amplatz Right		AR 1		672-110-00	558-110-00
		AR 1 SH	2 Side Holes	672-111-00	
		AR 1 LBT	Long BRITE TIP®	672-110-0L	
		AR 2		672-112-00	558-112-00
		AR 2 SH	2 Side Holes	672-113-00	
		AR 2 LBT	Long BRITE TIP®	672-112-0L	
Extra Back-up Right		XBR 1		672-122-00	
		XBR 1 SH	2 Side Holes	672-123-00	
		XBR 2		672-124-00	
Extra Back-up RCA		XB RCA		672-126-00	558-126-00
		XB RCA SH	2 Side Holes	672-127-00	
3 DRC— Williams Technique		3 DRC		672-130-00	558-130-00
		3 DRC SH	2 Side Holes	672-131-00	
NR—Noto Technique		NR 4		672-120-00	
Dedicated Transradial Shapes		Shape	Description	PRODUCT CODE	
Jean Fajadet Left		JFL		672-162-00	558-162-00
		JFL ST	Short tip	672-164-00	
Jean Fajadet Right		JFR		672-168-00	
		JFR ST	Short tip	672-170-00	
Radial-brachial Kimny Technique		RB		672-172-00	
Radial Bi-lateral		RBL 4		672-717-00	
		RBL 4.5		672-729-00	
		RBL 5		672-718-00	
Barbeau		BARBEAU		672-174-00	558-174-00
Right Radial		RRAD5			558-320-00
		RRAD6		672-300-00	
Right Radial Extra Back-up		RRAD6 XB		672-301-00	

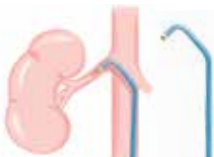
Multipurpose	Shape	Description	PRODUCT CODE	
			6F	5F
	MPA 1		672-270-00	558-270-00
	MPA 1 (125 cm)		672-270-25	
	MPA 1 SH	2 Side Holes	672-271-00	
	MPA 1 LBT	Long BRITE TIP®	672-270-0L	
	MPB 1		672-272-00	

Hockey Stick	Shape	Description	PRODUCT CODE				
			6F	5F			
				H-STICK		672-278-00	558-278-00
				H-STICK (90 cm)		672-278-90	
				H-STICK SH	2 Side Holes	672-279-00	
H-STICK LBT	Long BRITE TIP®	672-278-0L					

Left Coronary Bypass	Shape	Description	PRODUCT CODE				
			6F	5F			
				LCB		672-180-00	558-180-00
				LCB (90 cm)		672-180-90	
LCB SH	2 Side Holes	672-181-00					

Right Coronary Bypass	Shape	Description	PRODUCT CODE				
			6F	5F			
				RCB		672-182-00	558-182-00
				RCB (90 cm)		672-182-90	
RCB SH	2 Side Holes	672-183-00					

Internal Mammary	Shape	Description	PRODUCT CODE				
			6F	5F			
				IM		672-190-00	558-190-00
				IM (90 cm)		672-190-90	
IM SH	2 Side Holes	672-191-00					

Renal Curve	Shape	Description	PRODUCT CODE	
			6F	5F
	RDC		672-212-55	

ADROIT™ Guiding Catheter from Cordis

Speak with your Cordis representative about getting more out of your guiding catheter.

For healthcare professionals only. As part of the Cordis policy of continuous product development, we reserve the right to change product specifications without prior notification. Contact your Cordis sales representative for availability and ordering.
Important information: Prior to use, refer to the Instructions for Use supplied with this device for indications, contraindications, side effects, suggested procedure, warnings, and precautions.
 © Johnson & Johnson Medical NV/SA 2015 – EU736 05/15 - 2E-800-2074-4

